


waterpik

Ultra Cordless  
DENTAL WATER JET

**Quick  
Start  
GUIDE**

**1** Charge Unit for 16 Hours  
Prior to First Use.

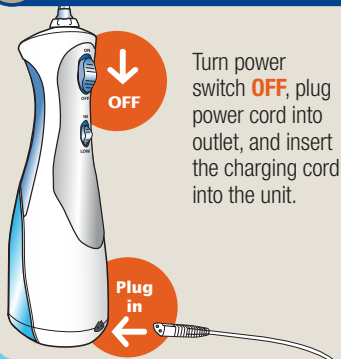


**2** Insert Tip:



**Caution:** Do not press tip eject button while unit is running.

**3** Fill the Reservoir:

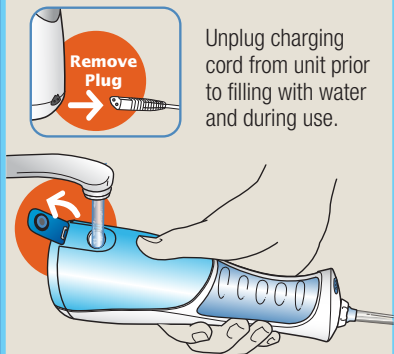

Lift the flip top on reservoir and fill with lukewarm water. Refill as needed.

**4** Adjust Water Pressure:





To increase water pressure move to **HI** setting.

*continued on back*

## 5 On/Off Procedure:

Always place tip in mouth before turning unit on.


Lean over sink and close lips enough to prevent splashing, while still allowing water to flow from mouth into the sink. Turn unit **ON**.

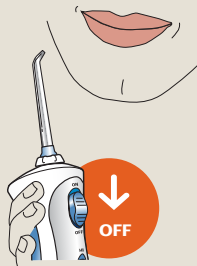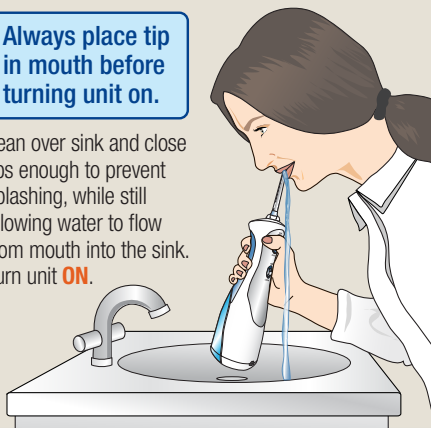


Always turn unit **OFF** before removing tip from mouth.

## 6 Technique Tips for Best Results:


For more details refer to the instruction booklet.

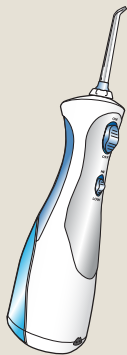

Keep unit upright during use to ensure continuous water flow.

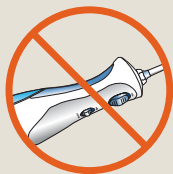

Aim the tip just above the gum line at a 90 degree angle.

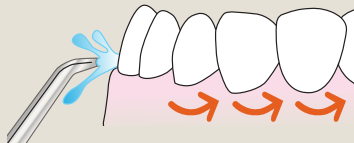

Starting with the back teeth, follow the gum line and pause briefly between teeth.

Continue until all areas around and between teeth have been cleaned.